


Product description:

A software for communication between a computer in PC class and standard Radwag weighing instruments. Main subject of the software is collecting the information on mass indicated by an instrument. The communication with the computer is realized through RS 232 interface.

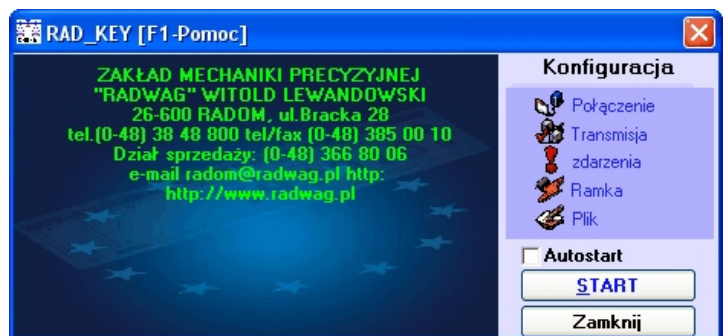

Designed for scales:

- with standard communication frame

Functions:

- collecting data from measurements and saving them into a file
- different ways of initializing the process of collecting data from an instrument to a computer (by a button in an instrument or computer keyboard)
- reading of transmission bits by RS 232 interface into a computer and placing the bits in place of cursor in the computer.

Computer requirements:

- processor 1,2 GHz
- free 500 MB on hard drive
- RAM 256 MB (recommended 512 MB)
- operation system Windows 2000/XP

Language versions:

Polish, Czech, German, English, French.

* NOTE! RADWAG does not sell software "in box". Boxes presented on our Web-site are of virtual characteristic only, and they are not accessible with purchased software. The customer receives purchased software on a CD in special protective cover and the user manual as a brochure.