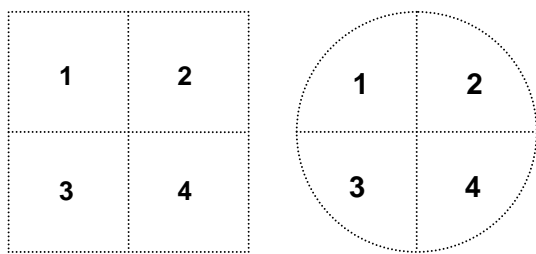


Examples of calibration of electronic weighing instrument with Max = 2100 g and d = 0,01 g

A. ERROR OF INDICATION

Load on the weighing pan	Error of indication	Extended uncertainty
500 mg	± 0 g	0,01 g
1 g	± 0 g	0,01 g
100 g	± 0 g	0,01 g
1000 g	± 0 g	0,01 g
2100 g	+ 0,01 g	0,01 g

B. ECCENTRICITY ERROR

Placing weight on pan	Load on the pan	Pole no.	Error of indication	Extended uncertainty
	700 g	1	+ 0 g	0,01 g
		2	+ 0 g	
		3	± 0 g	
		4	± 0 g	

C. DISPERSION OF READINGS

~50 % Max = 2000 g			~100 % Max = 4000 g		
Balance indication	Dispersion of indications	Extended uncertainty	Balance indication	Dispersion of indications	Extended uncertainty
1000,01 g	0,01 g	0,01 g	2000,01 g	0,01 g	0,01 g
1000,01 g			2000,01 g		
1000,00 g			2000,00 g		
1000,00 g			2000,01 g		
1000,00 g			2000,00 g		
1000,00 g			2000,00 g		
1000,00 g			2000,00 g		
1000,00 g			2000,01 g		
1000,00 g			2000,00 g		
1000,00 g			2000,01 g		
1000,01 g			2000,00g		